

The new Body Armor® Embrace Ankle Brace by DARCO.


Perfect stability along with controlled mobility!



- > post-traumatic and preventive care
- > anatomically shaped stirrups ensures a perfect fit
- > customizable
- > physiological correct movement of the foot
- > prevents the front shift of the talus



DARCO

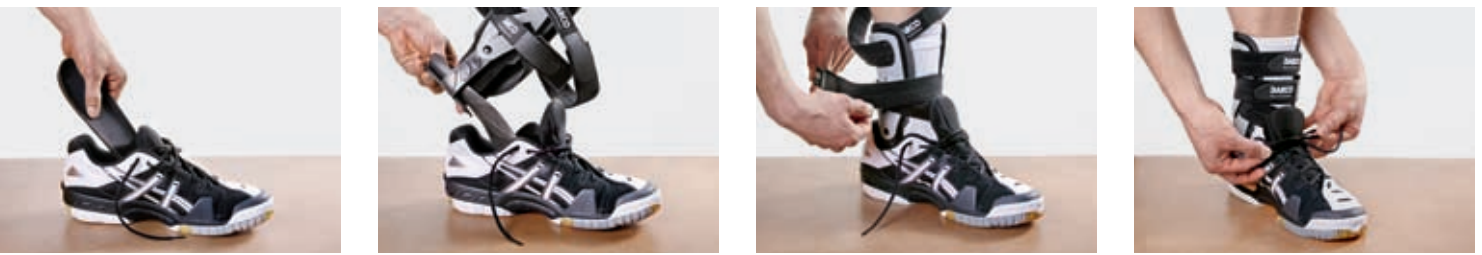

1. Elements of stability

Perfect fit with anatomically correct shaped elements of stability.

2. Orthosis sole

Additional space in the shoe due to shortened orthosis sole.

3. DARCO Y-closing-system

Individual strap customization with DARCO Y-closing-system.

4. Talus Shift Strap

Protection of the anterior talo-fibular ligament by means of the the Talus Shift Strap from plantar / lateral to dorsal / medial. (Prevents the front shift of the talus.)

5. Hinge joints

Physiological correct movement of the foot by parallel and thus anatomically exact positioning of the hinge joints.

Free dorsal extension and plantar flexion are possible.

Body Armor® Embrace

The Body Armor® Embrace is an ankle joint orthosis which is applied during the acute phase of injury as well as during the process of rehabilitation. The lower part of the insole (soft-pad) is mainly intended for padding the injured ankle and may be cut off at the marked segments after the ankle has healed. Two velcro fasteners warrant a smooth and easy but firm position of the foot. The Talus Shift Strap helps to prevent talar shift. The curved shape of the stirrup ensures a most comfortable fit and thus to an increased compliance of the user.

Indications

- > Distorsion trauma of the upper ankle joint
- > Insufficiency of the capsules and ligaments of the upper ankle joint
- > Prevention of chronic instability

Item-No.	Colour	Size	D / EU	UK
----------	--------	------	--------	----

left:

EB1B-L	Black-Grey	S	35-39	2,5-5,5
EB2B-L	Black-Grey	M	40-44	6,0-9,5
EB3B-L	Black-Grey	L	45-49	10-13

right:

EB1B-R	Black-Grey	S	35-39	2,5-5,5
EB2B-R	Black-Grey	M	40-44	6,0-9,5
EB3B-R	Black-Grey	L	45-49	10-13