

The original PegAssist™ System by DARCO for effective off-loading.



- > removable pegs for effective off-loading
- > shock-absorbing effect
- > available for the shoe models shown on reverse side
- > right and left fit
- > variety for OrthoWedge Light



DARCO


Pegs can be removed from specific areas of the insole in order to redistribute pressure. The cardboard cover then stabilizes the off-loading area.



The PegAssist™ System is available for all products listed here.

The 18 mm thick insole is layered with Poron®, Plastizote® and EVA to provide the perfect combination of shock absorption and moldability.

After removing the pegs from specific areas of the insole, the adhesive-backed cardboard cover is then applied to the complete surface of the insole. The cardboard cover acts to stabilize the remaining pegs.

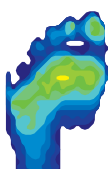

Standard Surgical Shoe



OrthoWedge



+ PegAssist™

MedSurg™

Item-No.	Size (D/EU)	
Women:		
PAS-W1	Small	34.0 – 36.5
PAS-W2	Medium	37.0 – 38.5
PAS-W3	Large	39.0 – 41.0
Men:		
PAS-M1	Small	39.0 – 41.0
PAS-M2	Medium	41.5 – 43.0
PAS-M3	Large	43.5 – 45.0
PAS-M4	X-Large	45.5 – 47.0

OrthoWedge Light (Variety)

Item-No.	Size (D/EU)	
PAS-STL0	X-Small	33.5 – 37.0
PAS-STL1	Small	37.5 – 39.5
PAS-STL2	Medium	40.0 – 41.5
PAS-STL3	Large	42.0 – 44.0
PAS-STL4	X-Large	45.5 – 47.0

All other models

PAS-0	X-Small	33.5 – 37.0
PAS-1	Small	37.5 – 39.5
PAS-2	Medium	40.0 – 41.5
PAS-3	Large	42.0 – 44.0
PAS-4	X-Large	45.5 – 47.0

