

Body Armor[®] Heel Reliever


Heel Protector

For
bedridden
patients



**Effective
pressure relief**
in the heel area.

The benefits

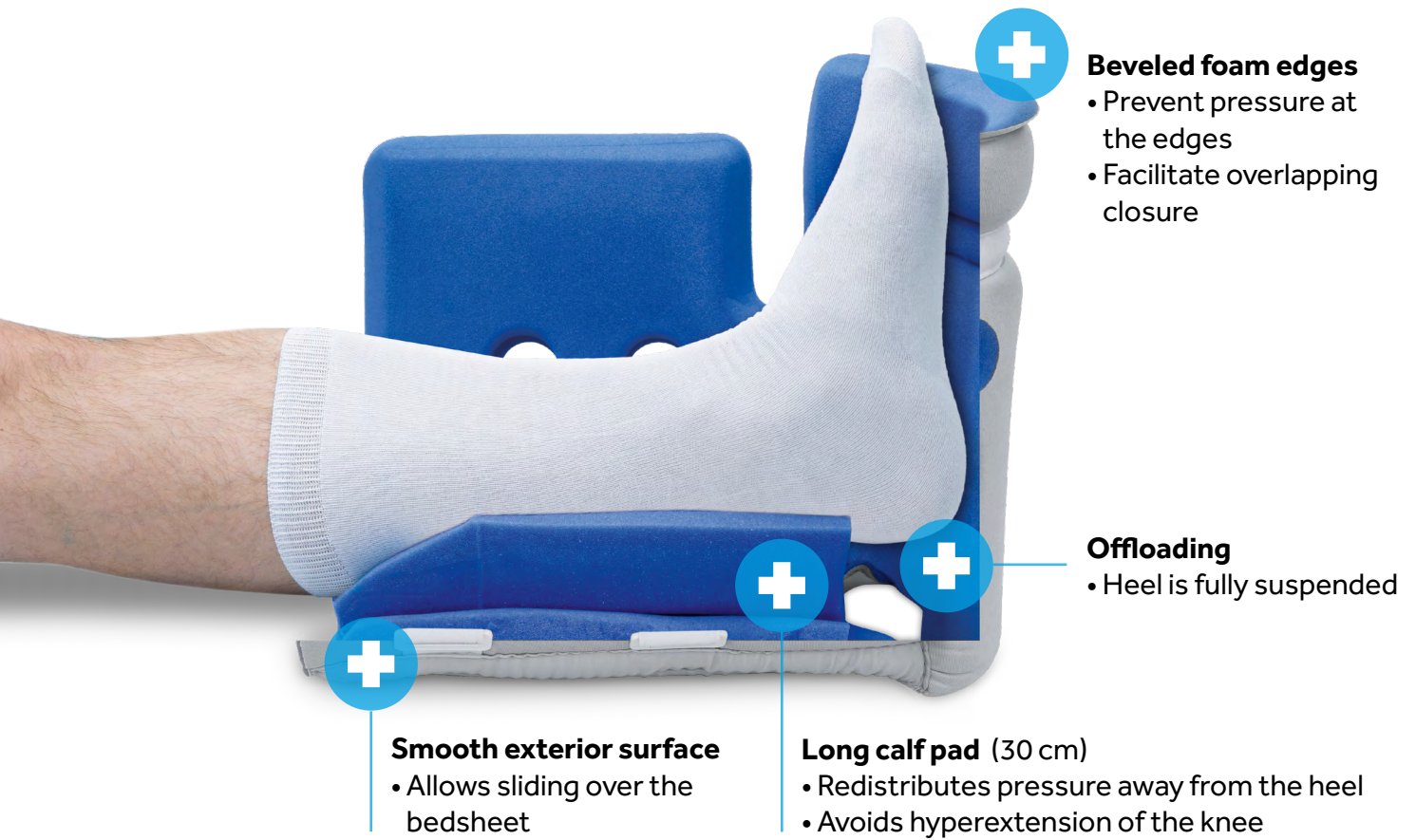

- Beveled foam edges**
- Prevent pressure at the edges
 - Facilitate overlapping closure

- Offloading**
- Heel is fully suspended

- Smooth exterior surface**
- Allows sliding over the bedsheet

- Long calf pad (30 cm)**
- Redistributes pressure away from the heel
 - Avoids hyperextension of the knee



- Pressure at the heel is redistributed to the Achilles tendon and the gastrosoleus complex of the calf
- The heel area is completely relieved of all contact pressure
- The heel is fully offloaded



- One-Size-Fits-All
- Adjusts to slim or larger calves
- The asymmetrical design allows overlapping of the side panels




- For visual inspection of correct heel offloading




- Suitable for use during negative pressure wound therapy (NPWT)

Indications


- For the prevention and healing of heel pressure ulcers as well as diabetic foot syndrome and decubitus
- Braden scale scores of 14 or less
- Diabetes with neurological sensory deficits
- Pedal pulse not palpable
- Paraplegia
- Prolonged immobility or coma

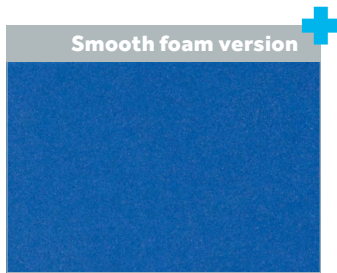

- Specifically designed for sensitive skin
- Particularly suitable for patients with diabetes
- Helps counteract edema formation

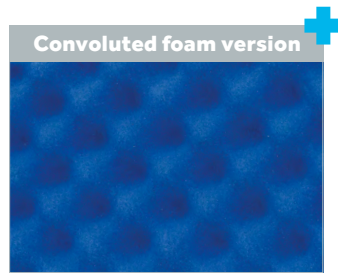


- Improves air circulation, providing a pleasant cooling effect

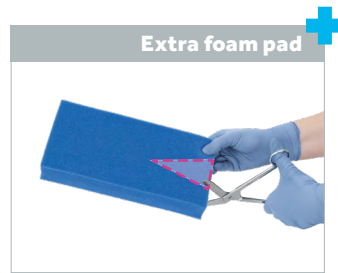

- Easy to cut to size
- Allows customized, pressure-redistributing positioning of the device



- Accommodates larger feet
- Maintains an orthograde position
- Prevents foot drop
- Prevents foot rotation
- Allows room of movement
- Helps prevent plantar flexion contracture

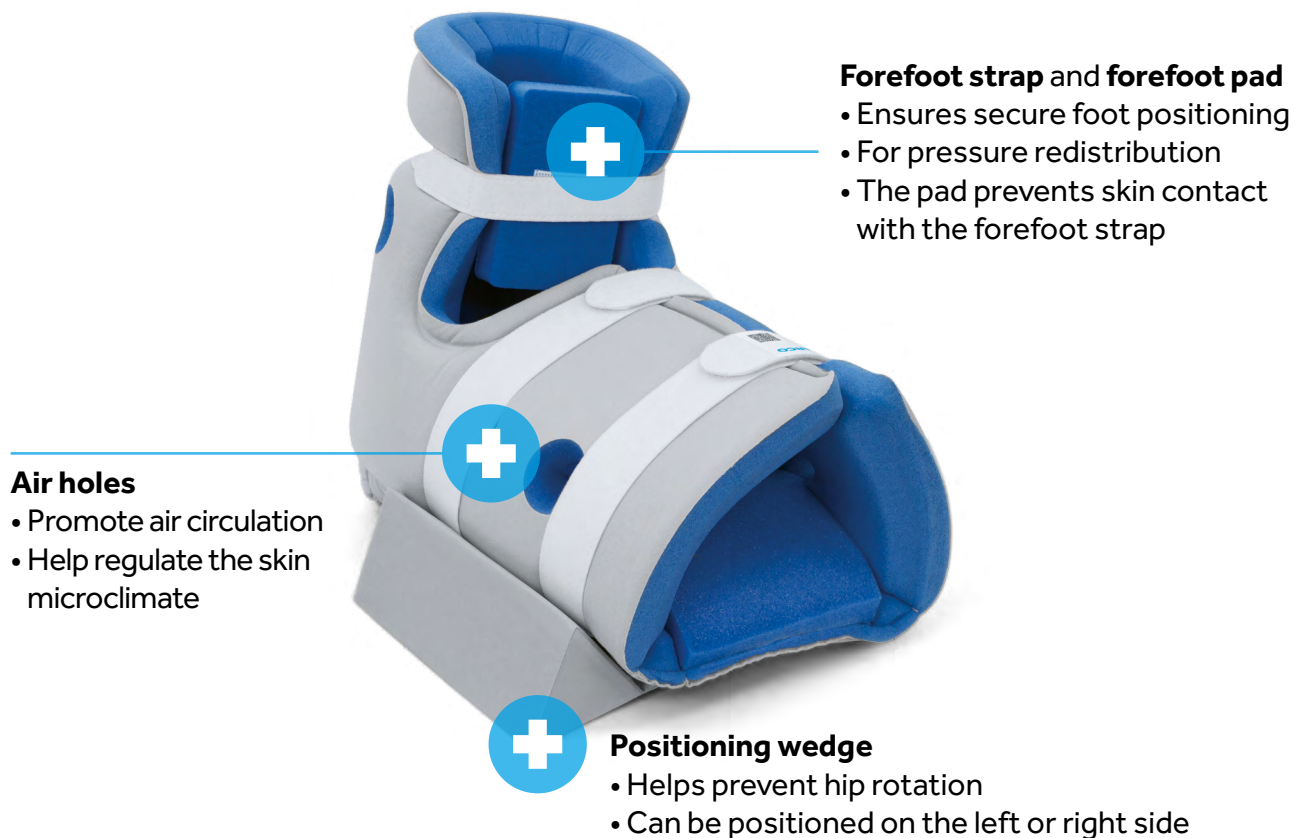







Sizes + measurements

Size	Shoe size	Circumference	Art. No.	Art. No.
			Convoluted	Smooth
Universal	up to 47	25 - 54 cm	BAHR-N	BAHR-F



Measure circumference

- Universal size 
- Right = Left 
- For diabetics 
- Low weight 
- Time saving 
- Washable at 60° 

Included in the package



Forefoot pad



Additional foam pad



Positioning wedge



Laundry bag



Product video



DARCO (Europe) GmbH

Gewerbegebiet 18 | 82399 Raisting | Germany | Tel. +49 8807 9228-0 | info@darco.de | www.darco.de

