

Body Armor[®] Pro Term

Foot Stump Orthosis

Interim care after Chopart and Lisfranc amputations



Right = Left



For diabetics



Time saving



90°
position



Perfect fit thanks to ease of
inflation and deflation



DARCO

The quick solution for mobilisation following forefoot amputation

Using the Body Armor® Pro Term Foot Stump Orthosis during post-op recovery, changes to the limb can be compensated using the soft lining with integrated pump. The BAPT enables individual adjustment to the changing conditions of the stump. The ankle-spanning shaft in a 90° position stabilises the upper and lower ankle. The ventral support strap and the variable/secure closure system enables circular, flat compression. The wedge-shaped outsole and insole in dorsal extension also ensures optimal rearfoot comfort.



Wide opening for easy foot insertion



Practical & secure closure system




Wedge-shaped insole system for optimal rearfoot weight-bearing

FEATURES

- > Interim care after forefoot amputations
- > Wedge-shaped outsole and insole in dorsal extension
 - Rearfoot weight-bearing
- > Shell construction: Foot to lower leg in 90° position
- > Wide opening
 - Easy foot insertion
- > Lining with integrated pump
 - Customised fit to stump conditions
- > Ventral support strap
 - Circular even compression
- > Secure closing system
- > Universal fit right and left
- > Colour: Black-Blue

SIZES AND MEASUREMENTS

ART. NO.	SIZE	STUMP LENGTH
BAPT2	M	16 cm
BAPT3	L	19 cm

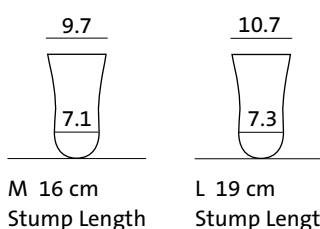

INDICATIONS

- > Postoperative care after Chopart / Lisfranc amputations
- > Diabetic foot syndrome
- > PAOD (peripheral arterial occlusive disease)
- > Trauma

COMPATIBLE WITH



WCS®
Wound Care Shoe