

MECRON ReadyDressment Strips

Elastic Stabilising Strips

Ready-to-use, elastic pre-cut tapes



Time saving



Postoperative tape stabilisation following:



DMMO



Claw / Hammer toe osteotomies



Chevron osteotomies

MECRON

Examples of use




Akin osteotomy

Surgical technique:

Osteotomy on the proximal phalanx of the big toe with a medial base wedge

Tape technique:

Tape running from the distal phalanx to the proximal phalanx on the medial side with gentle compression

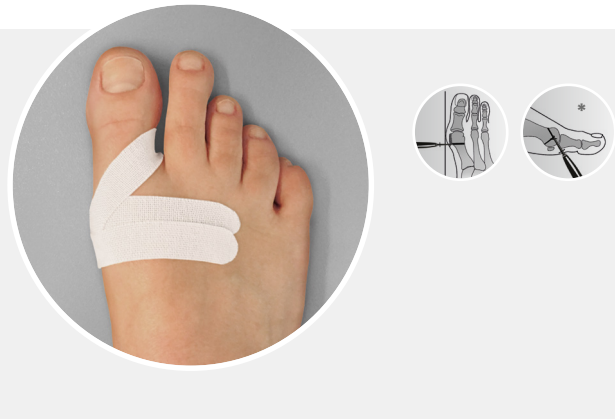

Chevron osteotomy

Surgical technique:

Osteotomy of the first metatarsal head to correct metatarsus primus varus in splayfoot deformities

Tape technique:

To secure the repositioning result with tape running from medial to lateral under slight tension

Indications

- > Postoperative and conservative use for the redressement and stabilisation of fingers and wrists, toes and ankles
- > Postoperatively following open foot procedures
- > Postoperatively for minimally invasive foot procedures:
 - Akin osteotomies
 - Chevron osteotomies
 - DMMO (distal metatarsal metaphyseal osteotomies)
 - Claw / Hammer toe osteotomies
 - Curly toe osteotomies
 - Digitus quintus varus osteotomies
 - Tailor's bunion osteotomies
 - Tenotomies

Features

- > Pre-cut in two different lengths
 - Individual anatomic fitting + control
 - Easy application + correction option = time saving
- > Provides a longer hold than wound closure strips
- > Early functional after-care possible
- > Elastic: Stretchable by approx. 60 %
 - No binding
- > Skin friendly dressing material
- > Latex-free
- > Can be worn for several days without restrictions
- > Dynamic stabilisation
- > Informative booklet with the first order (included)

Measurements

SIZE	WHITE ART. NO.	BLUE ART. NO.
Short 8 pieces, 10 × 1,4 cm (Length stretchable up to 16,5 cm)	MRSW-1	MRSB-1
Short BigPack 500 pieces, 10 × 1,4 cm	MRSW500-1	MRSB500-1
Long 4 pieces, 13,5 × 1,8 cm (Length stretchable up to 22 cm)	MRSW-2	MRSB-2
Long BigPack 500 pieces, 13,5 × 1,8 cm	MRSW500-2	MRSB500-2



Look into the booklet
(5-page preview | PDF):
"Application examples with the
MECRON ReadyDressment Strips"



* Acknowledgement: Photos from the book »Minimally Invasive Foot Surgery« courtesy of Dr. Mariano de Prado