

Movilo[®]

Therapeutic Comfort Shoe



3 new colours



So you can
enjoy walking
again!

Movilo® Fly

Berry



Black



Blue



Sand



Sky



Movilo® Street

Black



White



Olive



Pearl



Designed for more **comfort.**

The new comfort shoe from DARCO helps people who have problems with their feet or are generally on their feet a lot to get back to active, healthy and pain-free movement in everyday life and leisure. The Movilo® comfort shoe comes in two models: Fly and Street. Both are based on the same shoe last and there are only a few small differences between the two models. This way, we offer fashionable and comfortable shoes for people with foot problems - based on medical and therapeutic expertise, making life and walking easier.

Features



Outsole

- Back-shifted ball rocker
 - Forefoot off-loading
- Large abrasion-resistant rubber coating (3 mm)
 - Non-slip
- Shock-absorbing, force-absorbing sole construction
- Continuous and reinforced rigid shoe base
 - Forefoot off-loading
- Long and secure footprint
- Lower heel-to-toe drop (11 mm)



Insole System

- 2-part footbed with 8 mm total height
 - Anatomically moulded footbed (6 mm)
 - Moderate heel guide
 - Slight lengthwise and crosswise arch support
 - PU sole with smooth microfibre cover for good cushioning + easy foot insertion
 - Placeholder sole (2 mm)
 - For volume compensation
 - Additional shock absorption
- Replaceable
- Sufficient space (8 mm) for customised insoles
- Individually adjustable



Shoe Last Movilo®

- Foot-friendly last shape
- Ideal ball width
- Toe box height 21.5 – 27 mm (depending on size)
- Developed in cooperation with Wiest GmbH



Lacing

- Elasticated laces with quick-release fastener
 - Easy to put on and take off
- Non-elastic, flat laces (spare set enclosed)
 - Better support
- Additional lace hole
 - Secure fit



Shaft

- Air-permeable
- Reinforced toe
- Anatomical ankle cutout
 - Ideal rearfoot fit
- Optimally adapted heel cup
 - Substantial and cushioned heel support
- Loop on the heel
 - Easy to put on

Model Fly

- Loop on the shoe tongue
 - Easy to put on
- Without liner
 - Slim
 - Soft
 - Breathable


Model Street

- Stable upper material
- Soft liner
 - Good padding
 - Protection against pressure and external influences
 - Comfortable to wear

Indications

For the treatment of foot complaints:

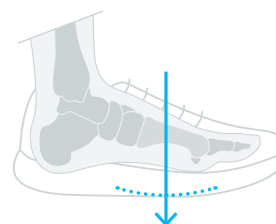


- Moderate foot / toe deformities such as hallux valgus, hammer and claw toes
- Hallux rigidus
- Metatarsalgia
- Morton's neuroma
- Rheumatoid arthritis
- Polyarthritis, arthrosis, gout
- Plantar fasciitis and heel spur
- Callus formation, panaritium prophylaxis
- Diabetes mellitus without PNP / PAOD (risk group 0 according to the German Diabetes Society)
- Diabetes mellitus without PNP / PAOD with foot deformity (risk group 1 according to the German Diabetes Society)
- Diabetes mellitus with sensory loss due to PNP / relevant PAD (risk group 2a according to the German Diabetes Society)



Back-shifted ball rocker





-  Exchangeable insole
-  For diabetics
-  Light weight

Sizes + Measurements

Movilo® Fly

EU Shoe size	Insole		Foot length	Ball girth	Toe box height	Art. No.				
	A	B				F	G	H	Berry ¹	Black ¹
36	24.0	9.0	22.5 - 23.1	22.6	2.15	MOFBYP36	MOFBAP36	MOFBUP36	MOFSDP36	MOFSKP36
37	24.6	9.1	23.1 - 23.7	23.0	2.20	MOFBYP37	MOFBAP37	MOFBUP37	MOFSDP37	MOFSKP37
38	25.3	9.2	23.8 - 24.4	23.4	2.25	MOFBYP38	MOFBAP38	MOFBUP38	MOFSDP38	MOFSKP38
39	26.0	9.3	24.5 - 25.1	23.8	2.30	MOFBYP39	MOFBAP39	MOFBUP39	MOFSDP39	MOFSKP39
40	26.6	9.4	25.1 - 25.7	24.2	2.35	MOFBYP40	MOFBAP40	MOFBUP40	MOFSDP40	MOFSKP40
41	27.3	9.6	25.8 - 26.4	24.6	2.40	MOFBYP41	MOFBAP41	MOFBUP41	MOFSDP41	MOFSKP41
42	28.0	9.8	26.5 - 27.1	25.0	2.45	MOFBYP42	MOFBAP42	MOFBUP42	MOFSDP42	MOFSKP42
43	28.6	10.1	27.1 - 27.7	25.4	2.50	MOFBYP43	MOFBAP43	MOFBUP43	MOFSDP43	MOFSKP43
44	29.3	10.2	27.8 - 28.4	25.8	2.55	-	MOFBAP44	MOFBUP44	MOFSDP44	MOFSKP44
45	30.0	10.3	28.5 - 29.1	26.2	2.60	-	MOFBAP45	MOFBUP45	MOFSDP45	MOFSKP45
46	30.6	10.4	29.1 - 29.7	26.6	2.65	-	MOFBAP46	MOFBUP46	MOFSDP46	MOFSKP46
47	31.3	10.5	29.8 - 30.4	27.0	2.70	-	MOFBAP47	MOFBUP47	MOFSDP47	MOFSKP47

Movilo® Street

EU Shoe size	Insole		Foot length	Ball girth	Toe box height	Art. No.			
	A	B				F	G	H	Black ¹
36	24.0	9.0	22.5 - 23.1	22.6	2.15	MOSBAP36	MOSWIP36	MOSOVP36	MOSPLP36
37	24.6	9.1	23.1 - 23.7	23.0	2.20	MOSBAP37	MOSWIP37	MOSOVP37	MOSPLP37
38	25.3	9.2	23.8 - 24.4	23.4	2.25	MOSBAP38	MOSWIP38	MOSOVP38	MOSPLP38
39	26.0	9.3	24.5 - 25.1	23.8	2.30	MOSBAP39	MOSWIP39	MOSOVP39	MOSPLP39
40	26.6	9.4	25.1 - 25.7	24.2	2.35	MOSBAP40	MOSWIP40	MOSOVP40	MOSPLP40
41	27.3	9.6	25.8 - 26.4	24.6	2.40	MOSBAP41	MOSWIP41	MOSOVP41	MOSPLP41
42	28.0	9.8	26.5 - 27.1	25.0	2.45	MOSBAP42	MOSWIP42	MOSOVP42	MOSPLP42
43	28.6	10.1	27.1 - 27.7	25.4	2.50	MOSBAP43	MOSWIP43	MOSOVP43	MOSPLP43
44	29.3	10.2	27.8 - 28.4	25.8	2.55	MOSBAP44	MOSWIP44	MOSOVP44	MOSPLP44
45	30.0	10.3	28.5 - 29.1	26.2	2.60	MOSBAP45	MOSWIP45	MOSOVP45	MOSPLP45
46	30.6	10.4	29.1 - 29.7	26.6	2.65	MOSBAP46	MOSWIP46	MOSOVP46	MOSPLP46
47	31.3	10.5	29.8 - 30.4	27.0	2.70	MOSBAP47	MOSWIP47	MOSOVP47	MOSPLP47

¹ Never Out of Stock (NOS): Permanently available

² Limited Edition (LE): Temporary available

Find the right size for you



All measurements are in cm and refer to the insole length. The indicated shoe sizes are for guidance only. The actual size depends primarily on the volume of the bandage or dressing.



DARCO (Europe) GmbH

Gewerbegebiet 18 | 82399 Raisting | Germany | Tel. +49 8807 9228-0 | info@darco.de | www.darco.de

